



1079 W / CO2 CRYO & EFC, Installation

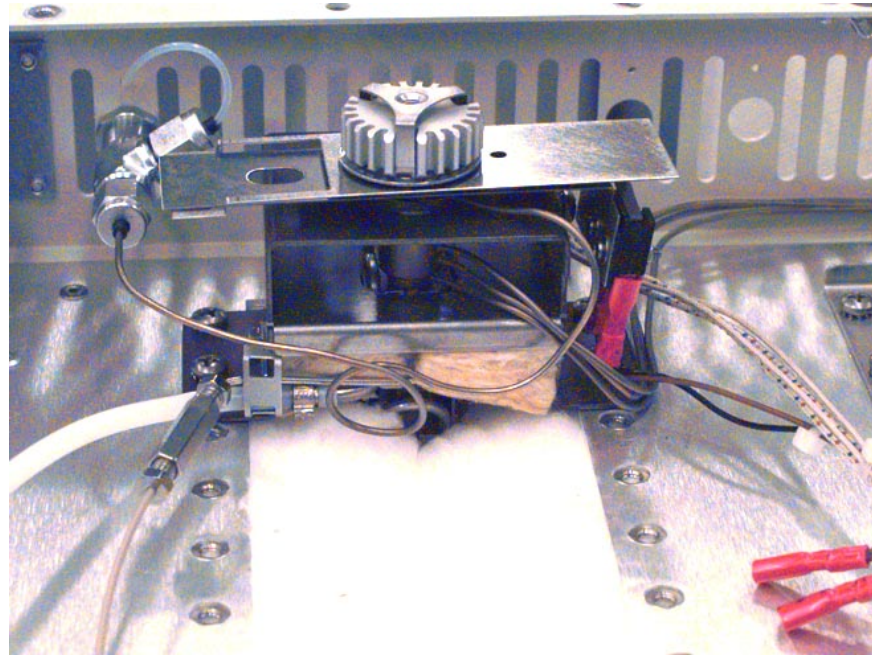

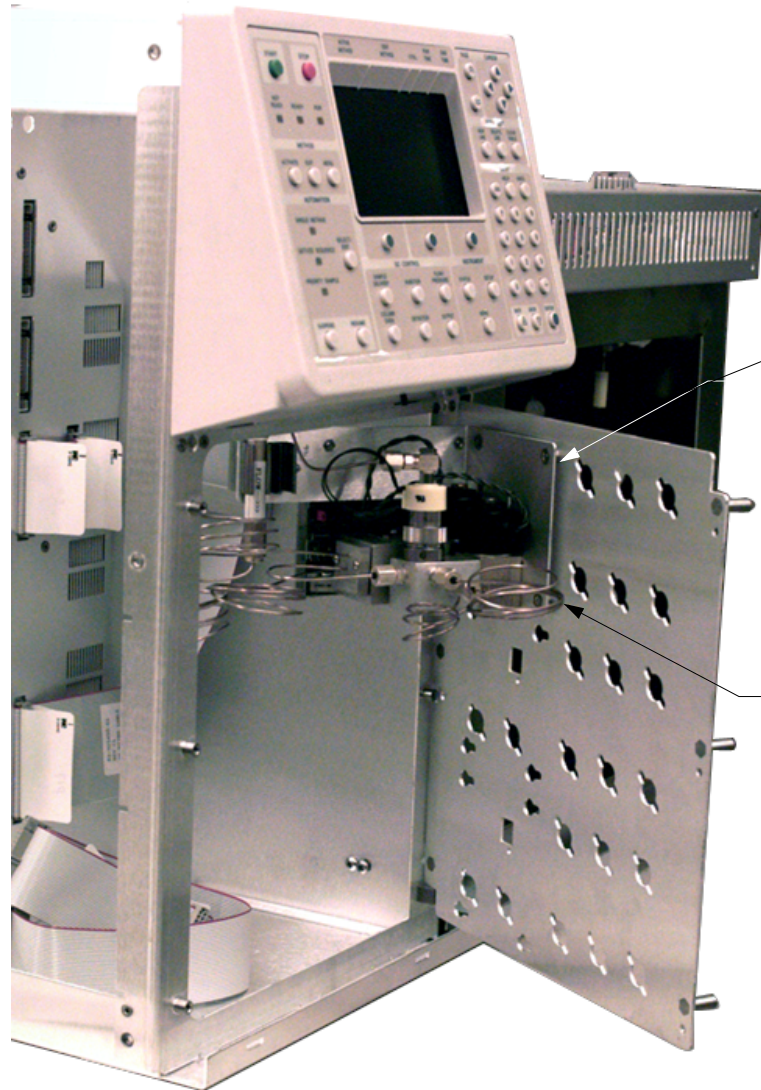

03-918532-01
VAR SEPTUM PURGE


12-901163-00
SCREW, 8-32 X 3/8"
EXT. LW 4EA. 8 TOTAL

28-211457-00
1/16-1/16 VALCO UNION

88-48071-00
1/18" INSULATION TUBING

03-925110-01
EFC TYPE 1

03-918853-00
PEEK TUBING

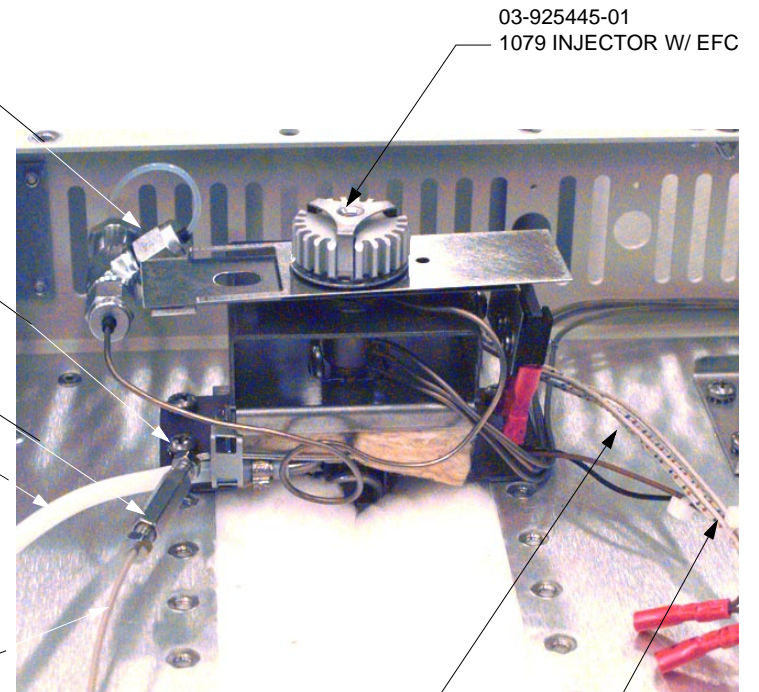

03-925445-01
1079 INJECTOR W/ EFC

03-925396-01
120V HTR. CART.

03-925396-02
220V HTR. CART.

03-925396-03
101V HTR.CART.

REF.

REF.


03-925374-01
TEMPERATURE PROBE

[CLICK HERE FOR
PROBE CONNECTION
TO J7 ON MAIN BOARD](#)

[CLICK HERE FOR
HEATER CONNECTION
TO J41 ON POWER BD.](#)

88-189609-00
TEFLON TAPE

28-694205-00
APADT, COMP 1/8" X 1/8"

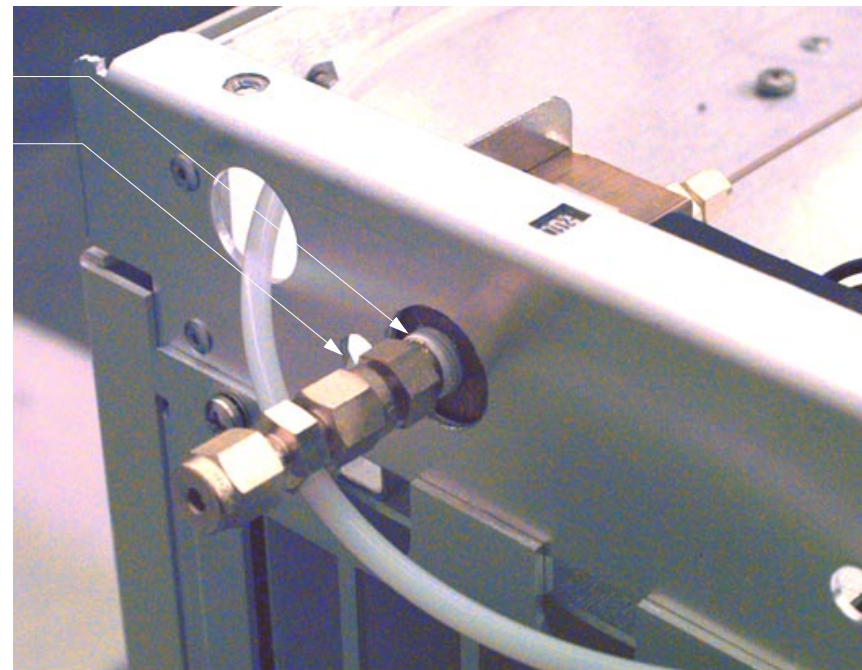


NOTE: PORT FACING INJECTOR
TO BE ARROW DIRECTION.

28-695110-00
ADAPT, COMP 1/16" X 1/8"

88-189609-00
TEFLON TAPE

03-925567-01
CO2 CRYO VALVE

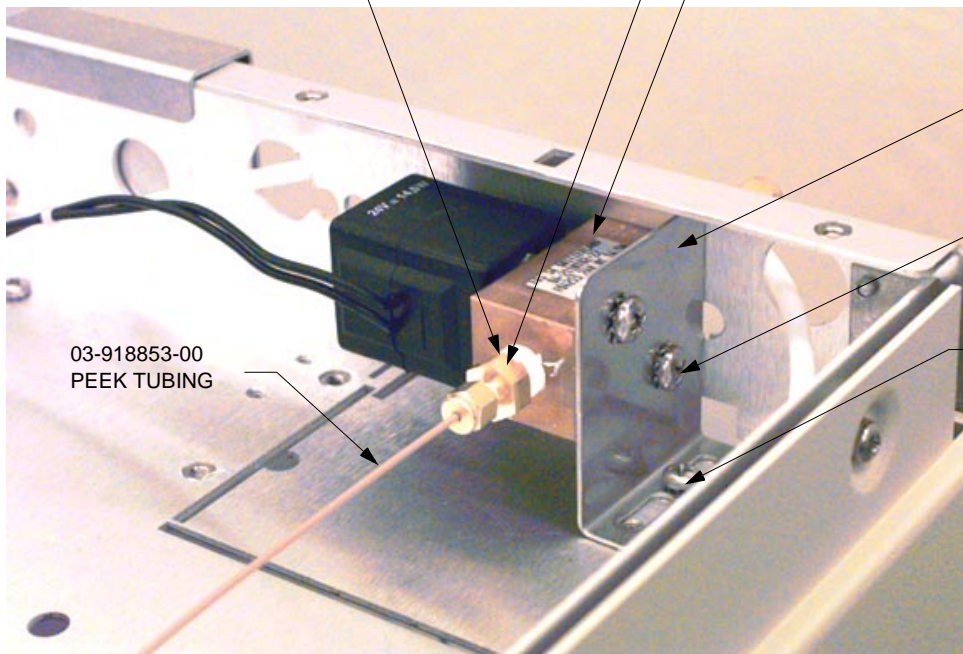


03-918853-00
PEEK TUBING

03-925562-01
CO2 CRYO BRKT

12-222060-04
SCREW 10-32 X 1/4" 2EA

14-312010-00
#10 EXT LOCK WASHER 2EA.

12-901163-00
8-32 X 3/8" EXT SIMS 2EA.

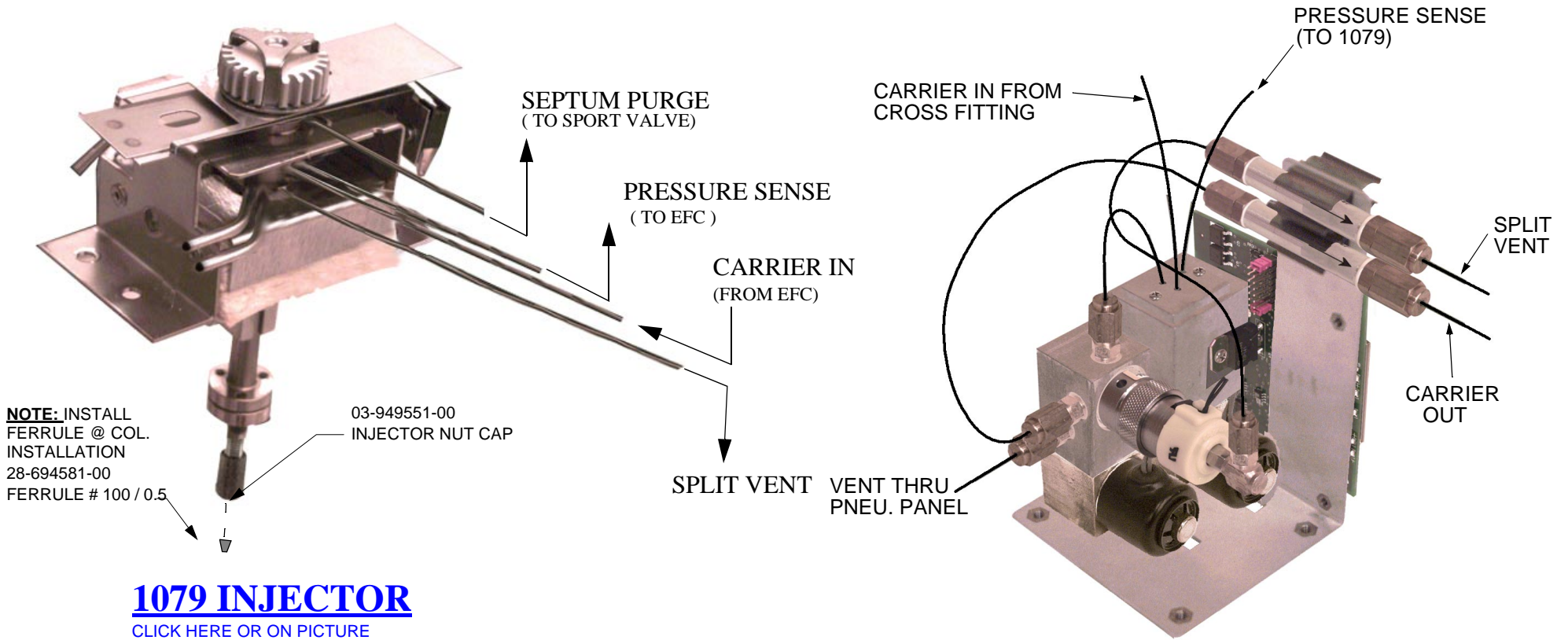

NOTE: INSTALL FERRULE @ COL. INSTALLATION 28-694581-00 FERRULE # 100 / 0.5

03-949551-00 INJECTOR NUT CAP

1079 INJECTOR

[CLICK HERE OR ON PICTURE TO SEE S/A EXPLOSION](#)

EFC TYPE 1

[CLICK HERE OR ON PICTURE TO SEE S/A EXPLOSION](#)

